

PureStation™ Dedicated Eye Sink

Cleaning Area & Cleaning Equipment All In One

The PureStation™ Dedicated Eye Sink serves as a designated cleaning area for channeled ocular devices. The PureStation™ Dedicated Eye Sink is permanently labeled with specific instructions for eye instruments, and can be used on a countertop, in a deep sink or on a mobile cart as an additional sink. The FlexiPump™ can flush up to five devices at once. A chemical mixing well and thermometer strip ensure that cleaning solutions are diluted and maintained at the proper temperature. Solution and rinse containers allow for the proper amount of fluid each time you pre-clean, and an accessory container keeps cleaning brushes and other accessories organized and secure. In addition, the sink is lightweight and easy to carry, with side handles that also serve as overflow prevention outlets.

Flush up to 5 devices at once

pure
PROCESSING

EYE Instruments Only

Discard H₂O used to clean or rinse after each use
Follow device manufacturer's instructions for use

www.pure-processing.com

(877) 718-6868



Caddy holds
parts and
brushes


Can be used on a
countertop or in a
mobile cart

Prevent loss of
small instruments
with drain catch

Properly mix rinse
solutions

Properly mix
cleaning solutions

Mixing well for
detergent

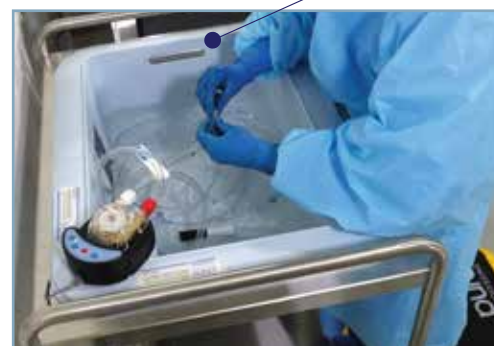

ANSI/AAMI ST79: 2017, M.2.2.2 Designated cleaning area and equipment A designated cleaning area and equipment dedicated to the cleaning of intraocular surgical instruments should be identified. Intraocular surgical instruments should be processed separately from general surgical instruments and equipment to reduce the potential for cross-contamination by material or residue from general surgical instruments.

pure
PROCESSING

WEB www.pure-processing.com

PHONE 877.718.6868

EMAIL sales@pure-processing.com