

# PureStation™ Dedicated Eye Sink

Cleaning Area & Cleaning Equipment All In One

The PureStation™ Dedicated Eye Sink serves as a designated cleaning area for channeled ocular devices. The PureStation Dedicated Eye Sink is permanently labeled with specific instructions for eye instruments, and can be used on a countertop, in a deep sink or on a mobile cart as an additional sink. The FlexiPump™ Independent Flushing System can flush up to five devices at once. Solution and rinse containers allow for the proper amount of fluid each time you manually clean, and an accessory container keeps cleaning brushes and other accessories organized and secure. In addition, the sink is lightweight and easy to carry, with side handles that also serve as overflow prevention outlets.

Flush up to 5 devices at once

**PURE**  
PROCESSING™

**EYE Instruments Only**

Discard H<sub>2</sub>O used to clean or rinse after each use. Follow device manufacturer's instructions for use.

[pure-processing.com](http://pure-processing.com)

(877) 718-6868



Water volume indicator and chemical mixing well help ensure accurate dilution rates


Caddy holds parts and brushes

Can be used on a countertop or in a mobile cart

Prevent loss of small instruments with drain catch

Cleaning solution caddy properly mixes cleaning solutions

Thermometer strip monitors temperature

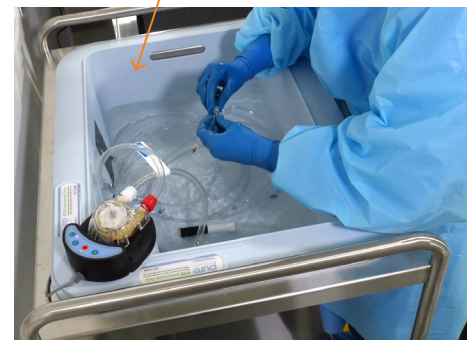

ANSI/AAMI ST79: 2017, M.2.2.2 Designated cleaning area and equipment A designated cleaning area and equipment dedicated to the cleaning of intraocular surgical instruments should be identified. Intraocular surgical instruments should be processed separately from general surgical instruments and equipment to reduce the potential for cross-contamination by material or residue from general surgical instruments.

**PURE**  
PROCESSING™

[pure-processing.com](http://pure-processing.com) | 877.718.6868