

New Generation of Pull Through Endoscope Brushes

Validated High-Level Cleaning Efficacy for Scratched Channels



Cleans endoscope channels in just
ONEPASS™ with validated performance!

SuperDART® three-stage design effectively eliminates debris and bioburden in just ONEPASS instead of 3!


- ◆ Leading Brush cleans debris from channel
- ◆ Squeegee conforms to channel design, pulling in enzymatic and removing bioburden
- ◆ Tail Brush removes remaining debris and bioburden from scratches and dents

*"Doesn't bring the
bio burden back
into the channel."*



ONEPASS does the work of three!

Scan to watch Product Video to learn
how SuperDART® works and why it's the
leading endoscope brush solution.

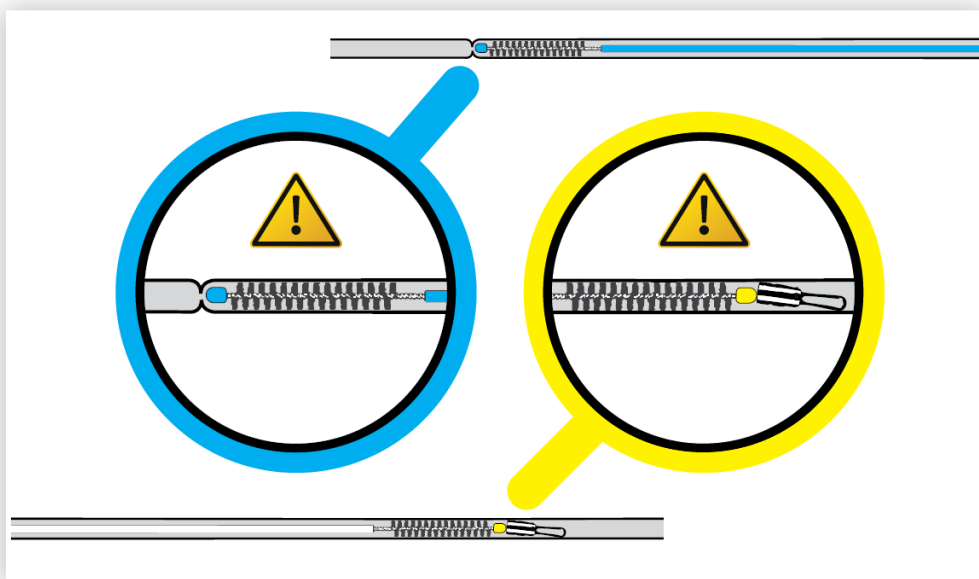

Color-coded
SmartSTOP™ prevents
accidental jamming and
detects channel flaws.

SuperDART® ONE PASS™ with SmartSTOP™

- ✓ Objective validated high-level of cleanliness
- ✓ Time savings with a single swipe
- ✓ Operators' fatigue reduction

Standard and Pull Through Brushes

- ✗ Subjective and uncertain cleanliness
- ✗ Multiple passes required
- ✗ Repetitive strenuous work

Provide your patients with the highest standard of care
by adopting the SuperDART® in your facility.

SuperDART® Part Number	Kit Part Number	Catheter Length (cm)	Catheter Diameter (mm)	Units/box
CHANNEL SIZE RANGE 2.0 - 2.8MM Ideal for Cystoscopes, Bronchoscopes, Intubation Scopes...				
9293	9399	240	1.6	25
9253	9355	122	1.6	25
9233	9333	85	1.6	25
CHANNEL SIZE RANGE 2.8 - 4.2MM Ideal for Colonoscopes, Gastrosopes, Duodenoscopes...				
9292	9188	250	1.7	25

Only use for the specified channel sizes
Always introduce the Leading Brush into the channels first
Always use wet, both the channel & the SuperDART®

